

# BABY SCALE
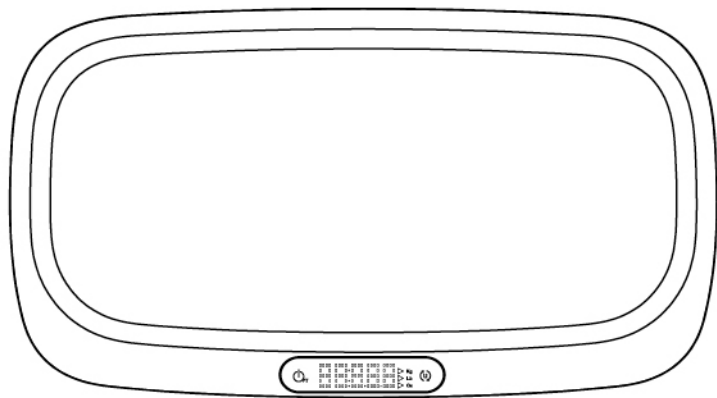
## USER MANUAL

# BABY SCALE

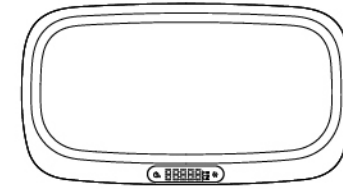Capacity  
28kg / 60lb

Thank you for purchasing the My Weigh® Baby Scale. Please read all operating instructions carefully before use. This electronic scale is a precision instrument. With normal care and proper treatment, it will provide years of reliable service. For more information please visit [www.myweigh.com](http://www.myweigh.com)

Never load the scale with more than the maximum capacity. Although the My Weigh® Baby Scale is designed to be extremely durable with extra overload protection built into the case, overloading will permanently damage it! Avoid any exposure to extreme heat or cold, your scale works better when operated at normal room temperature. Keep your scale in a clean environment. Dust, dirt, moisture, vibration, air currents and/or a close proximity to other electronic equipment can all cause an adverse effect on the reliability and accuracy of your scale. Handle with care. Gently apply all items to be weighed onto tray top. Avoid shaking, dropping or otherwise shocking the scale. Scales are delicate instruments and unlike cellular phones, scales have delicate sensors that determine how much an item weighs. If you drop or shock your scale, these sensors "feel" the shock and are sometimes destroyed. This happens with all digital scales. We design our scales to be as resistant to shock or drops as possible, however there is no way for us to protect 100% against load cell or sensor damage. Failure to follow these instructions will void your warranty.

Always allow the scale to acclimate to normal room temperature for at least one hour before use. Give your scale sufficient warm up time. Usually 30-60 seconds before calibration to give the internal components a chance to stabilize.

## PRECAUTIONS BEFORE USING THE BALANCE

1. Do not overload the scale in order to protect the sensor from damage.
2. Infants should always be placed on the center of the platform when being weighed.
3. Do not drop or shock the scale.
4. Clean the scale with a damp cloth, DO NOT immerse the scale in water or use chemical/abrasive cleaning agents.
5. Do not use propellant, abrasive or other chemicals to wipe the dirt in avoidance of discolour or malfunction.

## BATTERIES

Low Batteries & bad battery connections are the #1 cause of scale malfunction and inaccuracy! We test all of our scale returns from consumers and 60% of them are battery related problems. This sounds silly but it's true! A scale will perform poorly when it has low batteries. Use good quality batteries & replace them often (Remove the batteries if you plan to store the scale for longer than 14 days). If your scale simply won't turn on while on battery power, it is often caused by loose battery connections. Battery prongs (terminals) are made of metal and they have to be in contact with the batteries. You can use a paperclip to slightly bend the battery prongs to have a better connection. Some poorly designed batteries have recessed or partially obstructed battery terminals. This may cause your prongs to be touching the plastic housing of the battery instead of the metal of the battery terminal.

## Battery installation

- Press and lift open the battery cover located at the bottom of the unit.
- Insert 4 x AA batteries and make sure the polarity is correct (+) and (-).
- Close the battery cover until it clicks shut and the digits "0.00kg" are shown on the LED.

Note : If the battery symbol appears in the display, it means low battery power. It is time to replace the battery.

## OPERATION INSTRUCTIONS

### Select Measuring Unit

With batteries correctly installed, press UNIT button [⊕] in the scale to select measurement unit. The default measurement unit is "kg". You may press UNIT button [⊕] to choose between kilogram, pound and ounce.

### Daily Measurement

Turn on the scale and place the infant upon the scale until the LED displays the weight and stops blinking.

### TARE

TARE function allows you to weigh different items consecutively with out removing the items from the scale.

- When the scale is off, press the ON/OFF/TARE button [⊖] to go to the weighing mode. You can start to weigh when "0.00" is shown on the LED.
- Place the first item you wish to weigh onto the scale, the LED will show the weight of the first item.
- Press the ON/OFF/TARE button [⊖] to reset the display to zero.
- Now add another item to weigh onto the scale. The LED will show the weight of the newly added item.
- If you remove the added item from the scale, the LED will again show zero. But if you remove all items from the scale, the LED will show "-----".
- When it shows the "-----", it will automatically clear to zero and go to the normal weighing mode.

### CALIBRATION

Calibration may be required when the scale is first set up for use, or if the scale is moved to a different altitude or new location. This is necessary because the weight of a mass in one location is not necessarily the same in another location. Also, with time and use, mechanical deviations can occur.

How to calibrate: \*\*you must have an accurate 20 kg weight or combination of weights in order to calibrate\*\*

- Power on the scale and wait until it show zero.
- Press and hold ON/OFF [⊖] and UNIT button [⊕] at the same time for 3 seconds, display will show CAL, then several random numbers.
- Press UNIT button [⊕] one time, the display will show 20.00kg and flashing; put on 20kg test weight on the platform, then after stable, 20.00kg will stop flashing.
- Remove the weight from the platform and turn off the scale. The calibration process has been completed.

## FUNCTIONS

### Power

Pressing the POWER icon [⊖] will turn off the scale.

### Setting up to edit Auto Off and Brightness

- Turn on the scale and wait until "0.00kg" is shown.
- Press the UNIT button [⊕] for three times and then long press it one more time until you enter the Auto Off settings mode.
- You can now select in between the following Auto Off setting by using the ON/OFF/TARE button [⊖] to toggle between e.g. A-0, A-1, A-2, A-3 and confirm your selection by pressing the UNIT button [⊕].
- You enter the Brightness settings mode. You can now select in between the following settings by using the ON/OFF/TARE button [⊖] to toggle between e.g. L-0, L-1, L-2, L-3 and confirm your settings by long pressing the ON/OFF/TARE button [⊖] and turn off the scale.

### Auto Off

- Press A-0 to switched off completely auto off
- Press A-1 to put on auto off for 60 seconds
- Press A-2 to put on auto off for 120 seconds
- Press A-3 to put on auto off for 180 seconds

### Brightness

- Press L-0 to brighten LED light to level 0
- Press L-1 to brighten LED light to level 1
- Press L-2 to brighten LED light to level 2
- Press L-3 to brighten LED light to level 3

## CONNECTING THE APP TO THE BABY SCALE

- Turn on Bluetooth service on your device.
- Turn on Location service on your device. (Android only)
- Open My Weigh Baby app and proceed to the 'Weigh' screen.
- Turn on scale.
- Connection will happen automatically.
- If the scale disconnects press 'Tap to Reconnect' to attempt to re-establish the connection.

## HISTORY GRAPH

This screen displays the stored data for the current user. Weight data can be displayed over different periods of time. To see data for an individual date press the point on the graph that you wish to display. This will take you to a screen that displays the pictures, notes and weight data for that weighing event. Pressing the 'Delete' button will remove the data from your users profile and return you to the graph. More detailed instructions can be found by clicking the 'Help' icon on the 'History Graph' screen.

## PROFILES

The My Weigh Baby app allows for up to 5 user profiles. This allows you to track weight data separately for multiple users. To add a profile press the 'Plus' icon and enter the new users information before saving. Profiles may be edited or deleted by pressing the 'Options' button on the profile card. Please note that all weight data will be lost if a profile is deleted.

## SHARING TO SOCIAL MEDIA

When saving a new data point or when reviewing a previously saved point there is an option to share to social media. Pressing the 'Share' icon will take you to a screen with the data to be shared. Here you will find several options to share to common social media applications.

## SETTINGS

### Preferred Units

This setting allows you to choose the default format that weights will be displayed in, either Metric or US Customary.

### Automatically Speak Weight

When enabled this setting causes the 'Weigh' screen to read the weight aloud when the weight is settled on the scale.

## SPECIFICATIONS

Capacity	28kg / 60lb	Units	g, lb:oz, kg, lb, oz
Auto-OFF	A-0 = Off, A-1 = 60s, A-2 = 120s, A-3 = 180s		
Scale dimension	631mm x 347mm x 81mm (24"x 13" x 3")		
Cradle dimension	631mm x 347mm x 81mm (24"x 13" x 3")		
Scale weight	1.86kg		
Operating temperature	Optimum 5-40°C (41-104°F)		
Power Source	4 x C Batteries		
Tare Range	Up to scales maximum capacity		

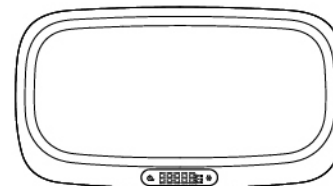
[www.myweigh.com](http://www.myweigh.com)

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